

# Direct Ophthalmic Instruments

Ophthalmoscope

Retinoscope

Otoscope

INSTRUCTIONS FOR USE



**Keeler**  
– A world without vision loss –



## CONTENTS

<b>1. INDICATIONS FOR USE .....</b>	<b>3</b>
<b>2. SAFETY.....</b>	<b>4</b>
2.1 PHOTOTOXICITY .....	4
2.2 WARNINGS AND CAUTIONS .....	4
2.3 CONTRAINDICATION.....	7
<b>3. CLEANING AND DISINFECTION INSTRUCTIONS.....</b>	<b>8</b>
3.1 STERILISATION .....	8
<b>4. INSTRUMENT HEADS.....</b>	<b>9</b>
4.1 OPHTHALMOSCOPES.....	9
4.2 LENS WHEEL.....	10
4.3 LENS RANGES .....	10
4.4 GRATICULE CONTROL .....	10
4.5 FILTER CONTROL .....	11
4.6 RETINOSCOPES .....	12
4.7 OTOSCOPES .....	12
4.8 BULB REPLACEMENT.....	14
<b>5. INSTRUMENT HANDLES .....</b>	<b>15</b>
5.1 HANDLE IDENTIFICATION .....	16
5.2 INSERTING/REPLACING BATTERIES .....	16
5.3 UPGRADE FROM BATTERY TO RECHARGEABLE HANDLES .....	16
5.4 BATTERY CONDITIONING .....	16
<b>6. GENMED WALL UNIT .....</b>	<b>17</b>
6.1 WALL MOUNTING.....	17
6.2 POWER SUPPLY ASSEMBLY .....	18
6.3 CONNECTING YOUR INSTRUMENT HEAD TO THE WALL UNIT HANDLE .....	18
6.4 DISPOSE-A-SPEC .....	19
<b>7. LITHIUM MINI CHARGER AND LITHIUM DOUBLE CHARGER.....</b>	<b>19</b>
7.1 POWER SUPPLY .....	19
<b>8. WARRANTY .....</b>	<b>20</b>
<b>9. SPECIFICATIONS AND ELECTRICAL RATINGS.....</b>	<b>20</b>
9.1 ELECTROMAGNETIC EMISSIONS.....	21
9.2 ELECTROMAGNETIC IMMUNITY .....	21
9.3 RECOMMENDED SAFE DISTANCES .....	23
<b>10. TECHNICAL SPECIFICATIONS.....</b>	<b>24</b>
<b>11. ACCESSORIES AND SPARES .....</b>	<b>25</b>
<b>12. PACKAGING AND DISPOSAL INFORMATION.....</b>	<b>27</b>

	Consult instructions for use		General warning sign
	Date of manufacture		Warning: Electricity
	Manufacturer's name and address		Warning: Floor level obstacle
	Country of manufacture		Warning: Non-ionizing radiation
	Waste Electrical and Electronic Equipment (WEEE) recycling		Warning: Optical radiation
	This way up		Warning: Hot surface
	Keep dry		Conformité Européene
	Fragile		Type B applied part
	Do not use if package is damaged		Class II equipment
	Temperature limit		Atmospheric pressure limitation
	Authorised representative in the European Community		Humidity limitation
	Catalogue number		Serial number
	Translation		Medical device

The Keeler Ophthalmic and Direct Instruments are designed and built in conformity with Directive 93/42/EEC, Regulation (EU) 2017/745 and ISO 13485 Medical Devices Quality Management Systems.

Classification: CE: Class I

FDA: Class II

The information contained within this manual must not be reproduced in whole or part without the manufacturer's prior written approval. As part of our policy for continued product development we the manufacturer reserve the right to make changes to specifications and other information contained in this document without prior notice.

This IFU is also available on the Keeler UK and Keeler USA websites.

Copyright © Keeler Limited 2021. Published in the UK 2021.

**Ophthalmoscopes:**

Pocket, Professional, Practitioner, Specialist, Standard

**Retinoscopes:**

Professional Combi, Spot, Streak

**Otoscopes:**

Fibre-Optic, Pocket, Professional, Standard

**Handles:**

Pocket, Slimline, GenMed Wall Unit

**Chargers:**

Lithium Duo Charger, Lithium Mini Charger

## 1. INDICATIONS FOR USE

These devices are intended to be used only by suitably trained and authorised healthcare professionals.



**CAUTION: Federal Law restricts this device to sale by or on the order of a physician or practitioner.**

### Intended use / purpose of instrument

The Keeler Ophthalmoscope is indicated for examining the posterior segment of the eye referred to as the fundus, in order to aid in the screening and diagnosis of retinal pathology including, but not limited to, diseases such as cataracts, papilledema, glaucomatous disc cupping, diabetic retinopathy, hypertensive retinopathy and retinal detachments. When set to a high power and magnification, it can also be used to examine the anterior segment of the eye, which includes the eyelids, cornea, sclera, conjunctiva, iris, aqueous, crystalline lens and anterior vitreous.

The Keeler Retinoscope is indicated for the objective assessment of refractive state of the eye. By observing the retinal red reflex, it is also a means to providing information regarding the visual system, for example media and lens opacities, significant ocular aberrations and accommodative status.

The Keeler Otoscope is indicated for the examination of the health of the external auditory canal, the tympanic membrane and the middle ear. Otoscopy can aid in the detection of ear conditions, including but not limited to earache, ear infection, hearing loss, ringing in the ears, inflammation and foreign bodies.

## 2. SAFETY

### 2.1 PHOTOTOXICITY



**CAUTION: The light emitted from this instrument is potentially hazardous. The longer the duration of exposure, the greater the risk of ocular damage. Exposure to light from this instrument when operated at maximum intensity will exceed the safety guideline after 4 hours 20 mins.**



**While no acute optical radiation hazards have been identified for Ophthalmoscope/ Retinoscope instruments, we recommend keeping the intensity of the light reaching the patient's retina to the minimum possible for the respective diagnosis. Children, people with aphakia and people suffering from eye conditions are most at risk. An increased risk may also occur if the retina is exposed to the same or a similar device with a visible light source within 24 hours. This applies, in particular, if the retina has been photographed with a flashbulb in advance.**

**Keeler Ltd shall on request, provide the user with a graph showing the relative spectral output of the instrument.**

### 2.2 WARNINGS AND CAUTIONS

Please note that the proper and safe functioning of our instruments is only guaranteed if both the instruments and their accessories are exclusively from Keeler Ltd. The use of other accessories may result in increased electromagnetic emissions or reduced electromagnetic immunity of the device and may lead to incorrect operation.

Observe the following precautions in order to ensure safe operation of the instruments.



#### WARNINGS

- Never use the instrument if visibly damaged and periodically inspect it for signs of damage or misuse.
- Check your Keeler product for signs of transport / storage damage prior to use.
- Do not use in the presence of flammable gases / liquids, or in an oxygen rich environment.
- US Federal Law restricts this device to sale by or on the order of a physician or practitioner.
- This device is intended to be used only by suitably trained and authorised healthcare professionals.
- This product should not be immersed in fluid.
- The power switch and mains plug are the means of isolating the device from the mains supply - ensure both the power switch and mains plug are accessible at all times.
- Do not position the equipment so that is difficult to press the power switch or remove the mains plug from the wall socket.



- Switch off the electrical supply and disconnect from the mains electrical supply before cleaning and inspection.

- If the product emits a strange odour, heat, or smoke, stop use immediately. The continued use of a damaged product or part may cause injuries.
- Do not touch the terminal contacts of the charging dock or hand unit, or the terminal contacts and the patient simultaneously.

**CAUTION**

- Use only genuine Keeler approved parts and accessories or device safety and performance may be compromised.
- Use only Keeler approved batteries, chargers, and power supplies as per the accessories listed in section 11.
- Backwards compatibility of the LED module has not been tested.
- The product has been designed to function safely when at an ambient temperature between +10°C and +35°C.
- Refraction stand variants or adaptors should only be used in combination with EN/IEC 60601-1 and EN/IEC 60601-1-2 compliant power supplies and devices.
- Keep out of the reach of children.
- To prevent condensation from forming, allow instrument to come to room temperature before use.
- For indoor use only (protect from moisture).
- There are no user serviceable parts inside. Contact authorised service representative for further information.
- Ensure device is securely held in docking station to minimise risk of injury or damage to equipment.
- Follow guidance on cleaning / routine maintenance to prevent personal injury / damage to equipment.
- Failure to carry out recommended routine maintenance as per the instructions in this IFU may reduce the operational lifetime of the product.
- At product end of life dispose of in accordance with local environmental guidelines (WEEE).
- To isolate the equipment, disconnect from the mains or switch off at the mains.
- The product and the ear specula are supplied non-sterile. Do not use on injured tissue.
- Use new or sanitised specula to limit the risk of cross-contamination.
- The disposal of used ear specula must occur in accordance with current medical practices or local regulations regarding the disposal of infectious, biological medical waste.

### Chargers



- Do not fit mains power adapter into a damaged mains outlet socket.



- Route power cords safely to eliminate risk of tripping or damage to user.
- Only Keeler handles with a red base can be used in the Keeler Lithium Chargers. Do not try to insert a Keeler handle with a blue base into the Keeler Lithium Chargers. Refer to Keeler handle and bulb identification.

### Direct Instruments

- When connecting instrument heads to handles please check that the voltage of the bulb in the instrument corresponds with the voltage of the handle.
- Care should be taken when fitting heads to handles not to trap skin between parts.
- Please ensure that the control is in the off position when the examination has been completed.
- Keeler Professional Retinoscopes contain strong magnets. Pacemakers and magnetically stored data will be affected or damaged by magnets.
- Strong magnetic fields may influence or distort sensitive electronic or mechanical test instruments. Very sensitive devices may even be destroyed. Always keep magnets at a safe distance from such devices.
- Do not use Keeler Retinoscopes or Ophthalmoscopes in ambient temperatures above 35°C.
- Disposable specula should not be used for insufflation testing.
- Plastic reusable specula will degrade if exposed to ultra- violet light, dry heat or gamma irradiation. These methods of sterilization must not be used.
- This device must only be used by clinicians trained in the use of ophthalmic devices.

### Batteries and LEDs

- Do not use a battery that is deformed, leaking, corroded or visually damaged. Handle a damaged or leaking battery with care. If you come into contact with electrolyte, wash exposed area with soap and water. If it contacts the eye, seek medical attention immediately.
- Ensure battery orientation is correct, or personal injury / damage to equipment may occur.
- Do not mix battery types.
- Do not attempt to charge non-rechargeable batteries.
- Do not charge battery in an environment where the temperature may exceed 35°C or fall below 10°C.
- When replacing rechargeable cell, turn handle off and insert new cell. Replace bottom cap, and place handle into charging well.
- If a short circuit occurs, reactivate the battery by placing the handle in the charger until the



LED flashes. This is a built-in protection device to protect the battery from damage.

- Dry cell batteries should be removed if your instrument is not to be used for long periods.
- Do not disassemble or modify the battery. There are no serviceable parts inside
- Do not dispose of battery in fire, puncture or short circuit.
- Dispose of batteries in line with local environmental regulations.
- Tape over battery contacts to avoid short circuiting during disposal.



- After removal of the battery do not touch the battery contacts and the patient simultaneously.



- Note: Lithium Ion and NiMH batteries contain no toxic heavy metals such as mercury, cadmium, or lead.



- Do not exceed maximum recommended exposure time.

- Always ensure that the handle rheostat is turned off before fitting an instrument head or changing a bulb.



- Bulbs / LED's can reach high temperatures in use – allow to cool before handling. The ophthalmoscope and retinoscope should not be continuously switched on for more than 15 minutes. If they are in the charging position or left on for 15 minutes or longer then they must be switched off and left to cool for at least 10 minutes before the next use.

- Care should be taken when handling halogen bulbs. Halogen bulbs can shatter if scratched or damaged.



- After removal of the bulb / LED do not touch the bulb / LED contacts and the patient simultaneously.

- Refer to the instructions on page 14 for bulb replacement.

### 2.3 CONTRAINDICATION

There is no restriction to patient population this device can be used with other than those outlined in the contraindications stated below.

Due to high illumination levels, the Ophthalmoscope and Retinoscope can cause some discomfort in photophobic patients.

Mydriatic agents when used in retinoscopy and ophthalmoscopy can cause temporary symptoms of blurred vision and photophobia. Adverse reactions to mydriatic drops are rare.

There are very few risks associated with otoscopy. Some patients may report a slight discomfort during the procedure, especially during the insertion of the speculum into a swollen and inflamed ear canal. If the plastic tip of the Otoscope is not replaced or cleaned properly, infection can spread from one ear to the other.

### 3. CLEANING AND DISINFECTION INSTRUCTIONS



**Before any cleaning of the instrument or the base unit, ensure the power lead is disconnected.**

Only manual non-immersion cleaning as described should be used for this instrument. Do not autoclave or immerse in cleaning fluids. Always disconnect power supply from source before cleaning.

1. Wipe the external surface with a clean absorbent, non-shedding cloth dampened with de-ionised water / detergent solution (2% detergent by volume) or water / isopropyl alcohol solution (70% IPA by volume). Avoid optical surfaces.
2. Ensure that excess solution does not enter the instrument. Use caution to ensure cloth is not saturated with solution.
3. Surfaces must be carefully hand-dried using a clean non-shedding cloth.
4. Safely dispose of used cleaning materials.

#### 3.1 STERILISATION

Plastic reusable Specula will degrade if exposed to ultra-violet light, dry heat, or gamma irradiation. These methods of sterilisation must not be used.



1. Re-useable specula should not be re-used if visibly contaminated with cerumen, ear drainage or blood. Dispose of safely.
2. Manually clean all surfaces of the units using a suitable brush and de-ionized water/ detergent solution (2% detergent by volume). Ensure that hinged versions of Specula are cleaned in both open and closed positions. Ensure all crevices are accessed. Solution can be heated to no more than 35°C.
3. Carefully examine to ensure that all visible contamination has been removed.
4. Safely dispose of used cleaning materials.
5. Sterilise using a validated steam steriliser complying with BS 3970 or equivalent standard. Operating cycle conditions as below: 134 - 138°C sterilising temperature at 2.25 bar operating pressure for minimum of 3 minutes hold time.



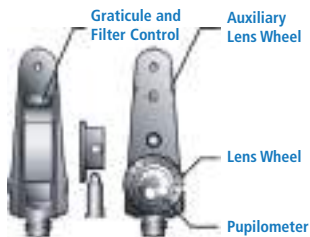
6. Following cleaning and/or sterilisation processes inspect the device to ensure all visible soil has been removed and the device operates as intended and is suitable for its intended use. Do not use if damaged. Dispose of safely.
7. The useful life of the device is determined by the wear and damage during use.

**Disposable Specula** – use once only and dispose of safely.


## 4. INSTRUMENT HEADS

## 4.1 OPHTHALMOSCOPIES


## Specialist




## Standard

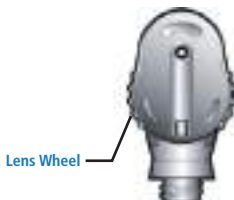

## Pocket

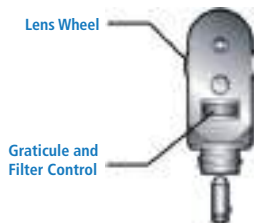

## Practitioner / Professional

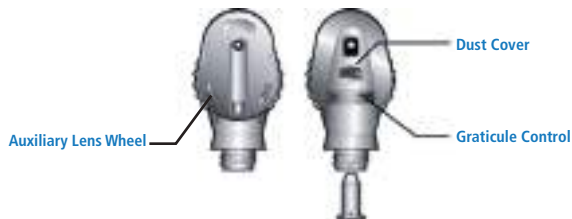

## 4.2 LENS WHEEL

The lens wheel is rotated to select the required lens. Lens powers are displayed in the viewing window as follows:

Black = (+) power lenses

Red = (-) power lenses

### Auxiliary lens wheel

Swings in +/- 20 in one Diopter step\* (\*Professional only).

### Specialist Auxiliary lens wheel

Rotate to align +10, +15, +30/ -10, -15, -30 Diopter lenses.

## 4.3 LENS RANGES

### Specialist

+44D to -45D in single Diopter steps

### Practitioner and Standard

+40D to -25D

### Professional

+29D to -30D in single Diopter steps

### Pocket

+20D to -20D

## 4.4 GRATICULE CONTROL

The graticule control is used to select the required beam for examination. The choice of graticules is as follows.



### Wide Angle

Illuminates the largest area of fundus for the best possible general diagnosis through a dilated pupil.



### Intermediate

Permits easier access through an un-dilated pupil in peripheral examination. Particularly useful in paediatric examination.



### Macular

Designed specifically for viewing the macular area of the fundus. Reduces pupillary reaction and improves patient comfort.



### Slit

Used primarily to determine retinal elevations and depressions but may also be used to assess anterior chamber depth.



### Glaucoma






























Projects a graticule onto the retina to assess the optic disc/cup ratio as an aid to glaucoma diagnosis and monitoring.



### Fixation Cross

Projects a graticule on to the retina for assessment of the degree and direction of eccentric fixation. This is particularly useful when examining children.

The Graticule Range for each ophthalmoscope is as follows:

Specialist							
Professional							
Practitioner							
Standard							
Pocket							

#### 4.5 FILTER CONTROL

The filter control\* is used to select the required filter.

(\*Professional/Practitioner/Standard only.)

#### Filter Applications



##### Red Free (Green filter)

Is used to examine the blood vessels in fine detail. The green filter blocks red rays showing blood vessels as black against a dark green background. This filter is particularly useful for diabetic retinopathy.



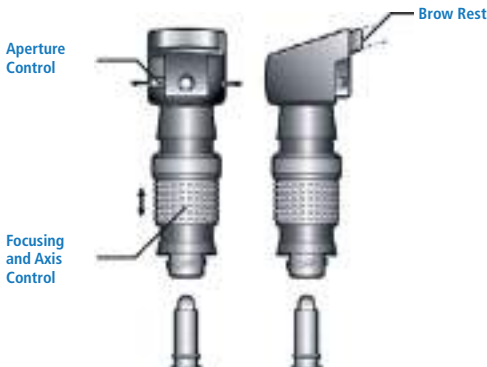
##### Cobalt Blue\*

Is used in conjunction with fluorescein dye for the detection and examination of corneal scars and abrasions. (\*Practitioner and Specialist only).

#### Pupillometer\*

Hold the pupillometer adjacent to the patients eye to estimate pupil size. 1=1mm. The range is 1mm to 8mm. (\*Applies to Specialist only.)

## 4.6 RETINOSCOPES



### Focusing and axis control (Streak)

The vergence is altered by sliding the focusing control up and down as indicated. In the top position the effect is a concave mirror. Mid position produces a streak behind the patient. The mid position is used to determine the presence and axis of any astigmatism. In the bottom position the effect is a divergent plane mirror effect. Refraction is normally performed between the mid position and the bottom position. The focusing and axis control can be rotated continuously in any direction.

### Focusing and axis control (Spot)

The vergence is altered by sliding the focusing control up and down as indicated. For all positions, the effect is a plane mirror effect.

### Brow Rest

The Keeler Retinoscope is supplied with a choice of brow rests to accommodate spectacle wearers. To interchange the brow rest, disconnect and attach as indicated.

### Aperture Control

The aperture control has two positions. To change from the large to the small aperture slide the control from left to right as indicated.

## 4.7 OTOSCOPES

Five permanent Specula are provided with each Otoscope / Set. The diameters are as follows: 2.5, 3.5, 4.5, 5.5 & 8mm. These are attached to the Otoscope head as shown in the following diagrams.

**Standard / Pocket****Fibre-Optic****Practitioner****Disposable Specula**

Disposable specula can be fitted to the Standard, Practitioner, Fibre-optic and Pocket Otosopes.

**Pneumatic Testing**

An insufflation tube can be fitted to your Otoloscope to enable you to carry out pneumatic testing.

For Practitioner, Standard, Pocket and Fibre-optic Otoloscopes, attach the insufflation adaptor into port. The insufflation tube can then be attached to this.

An Insufflation adaptor is also available for Practitioner as shown above.

**Minor Surgical Procedures**

Should you wish to use surgical instruments for minor procedures the following notes may be of assistance.

**Standard and Pocket Otoloscopes**

The magnifier can be removed to allow the introduction of surgical instruments.

**Fibre-Optic / Practitioner**

The Fibre-Optic magnifier can be moved to one side or removed completely to aid the introduction of surgical instruments.

## 4.8 BULB REPLACEMENT

Bulbs / LED's can reach high temperatures in use – allow to cool before handling.



- Always ensure that the handle rheostat is turned off before fitting an instrument head or changing a bulb.



- Care should be taken when handling halogen bulbs. Halogen bulbs can shatter if scratched or damaged.
- After removal of the bulb / LED do not touch the bulb / LED contacts and the patient simultaneously.
- Keeler bulbs can only be used in the instrument for which they are designed – refer to part number list in section 11. Ensure the replacement bulb is the correct voltage. See base of bulb.

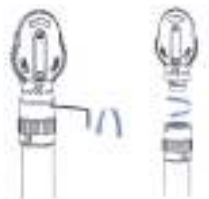
Blue = 2.8V for dry cell battery handles.

Red = 3.6V for rechargeable handles.

Black = LED.



- Loosen the set screw securing the instrument head to the handle. (GenMed Wall Unit only)
- Remove the head by holding it horizontally with one hand while rotating the handle counterclockwise with the other.
- Take care to ensure the battery / bulb does not drop out when the head and handle are separated.
- Remove the faulty bulb and dispose of in accordance with local environmental regulations.
- Replace the bulb with one of the correct voltage and type. Ensure that the location key is aligned with the aperture in the instrument head.
- Refit the handle to the head by rotating it clockwise while horizontal. If required, secure the head in place with the set screw provided. (GenMed Wall Unit only)





## 5. INSTRUMENT HANDLES

### Slimline



### Pocket



### Connection of the instrument heads to the handle

The connection between the instrument head to the handle is a screw thread. To connect our instrument head connect as shown and rotate in clockwise direction. Ensure the connection between the head and handle is positive.



### Compatibility

The Keeler Specialist, Professional, Standard and Practitioner Ophthalmoscopes and Keeler Retinoscopes are compatible with Keeler 2.8V and 3.6V Keeler handles.

The Keeler LED module is only compatible with Keeler 2.8V and 3.6V slimline handles.

### On / Off brightness control

To switch the instrument on, rotate the brightness control as indicated to the right.

To switch off the instrument, rotate the brightness control as indicated to the left.

Keeler Slimline Handles have a power indicator. This will show if the instrument is on or off.

Silver = off



Off

Red = on



Half On



On

## 5.1 HANDLE IDENTIFICATION

Keeler slimline handles are colour coded to allow you to distinguish between a dry cell battery handle (2.8v) and a rechargeable handle (3.6v).

The handles and Keeler bulbs are colour coded as follows:

Blue base = 2.8V for dry cell batteries.

Red base = 3.6V for rechargeable batteries.

Black base = LED for dry cell and rechargeable batteries.



- Please ensure when replacing batteries and bulbs that the voltage corresponds to the handle.

Disconnect from charger prior to removing instrument head.

Dispose of old batteries safely.

## 5.2 INSERTING/REPLACING BATTERIES

Unscrew battery cap, insert batteries, and replace battery cap as shown on page 15.



- Please note Keeler rechargeable handles are normally supplied complete with a rechargeable battery (3.6V).

### Dry cell batteries

The following dry cell batteries should be used:

- Keeler Pocket Handle – 2 x AA size dry cell batteries – Duracell MN 1500 or equivalent.

## 5.3 UPGRADE FROM BATTERY TO RECHARGEABLE HANDLES

Your Keeler 2.8V slimline (blue base) dry cell battery handle can be upgraded to a 3.6V (red base) rechargeable handle. Refer to section 11 for details of part numbers required.

Please note the bulb in your instrument will also need to be upgraded from 2.8V to 3.6V.

### Battery charging



- Do not attempt to charge non-rechargeable batteries.

## 5.4 BATTERY CONDITIONING

Your Keeler rechargeable batteries need to be conditioned to ensure you achieve the maximum life from the product. Follow the conditioning instructions as indicated.

### Step 1

Fully charge your new Keeler rechargeable battery. This will take approximately 15 hours.

### Step 2

Use the instrument without recharging until the battery is completely empty.

### Step 3

Once empty recharge the battery until full. This will take approximately 15 hours.

Repeat steps 1, 2 and 3 three times, i.e. fully charge and discharge the battery three times to complete the conditioning process. Once you have conditioned your batteries as described above you may place your instrument in the charger when not in use between examinations.

### Charger Compatibility



- Keeler Rechargeable Handles can be used in the following Keeler chargers only:
  - Keeler Mini charger
  - Keeler Duo charger

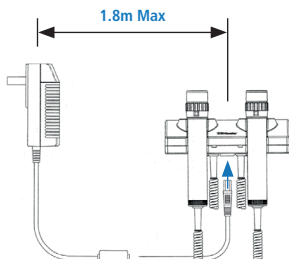


- Note: Handheld diagnostic instruments can become hot during use and charging.


## 6. GENMED WALL UNIT

### 6.1 WALL MOUNTING

Check the distance from the wall socket to the intended mounting position.





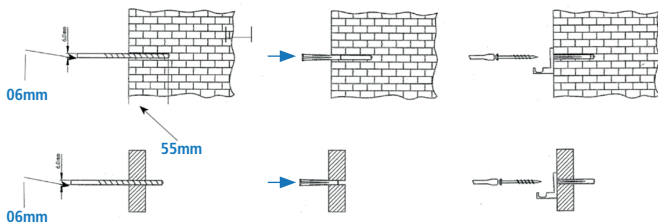For the Gen Med Wall Units drill two holes  $\text{Ø}6\text{mm}$  x 55mm deep and 100mm apart.



For the Dispenser Unit drill an additional two holes 249mm below the existing holes.

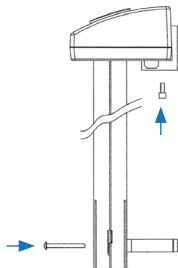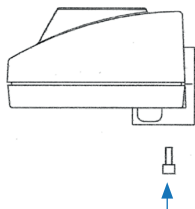Secure the GenMed Wall Unit and Dispenser Unit as shown.



## 6.2 POWER SUPPLY ASSEMBLY

### Set Plug

Replace the blanking plate with the appropriate mains plug adaptor if required, or use IEC 60320 TYPE 7 connector (not supplied).

**Please note:**



- This equipment may be affected by electromagnetic interference.
- Other electrical equipment in the close vicinity may also be affected by the GenMed Wall Unit.
- If such effects are suspected, switch off the offending equipment.

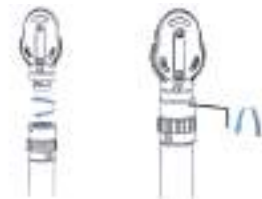
## 6.3 CONNECTING YOUR INSTRUMENT HEAD TO THE WALL UNIT HANDLE

The instrument head should be screwed positively onto the handle as shown.

As an additional security measure, instrument heads may be locked onto Keeler cord handles by tightening the built in screw with the hexagonal key provided.

To use the required instrument, remove the relevant handle from its cradle as shown.

A yellow light (LED) will illuminate when a cord handle is removed from its cradle. This will occur whether or not an instrument head is fitted.



When the instrument is no longer required always ensure that the handle is returned correctly to its cradle and that the LED goes out.

Only one handle can be used at a time. Replace the handle before using the other instrument.

Refer to the instructions in section 5 for information on the controls and operation of Ophthalmoscope, Otoscope and Retinoscope heads.

## 6.4 DISPOSE-A-SPEC

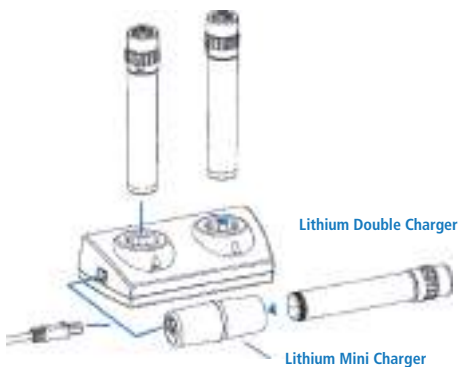
To dispense a speculum, simply grasp the end of the required speculum and gently pull vertically. When a dispenser tube is empty, reorder speculum using the order form EP59-48483.

Remove the lid from the unit and refill the required tube.

## 7. LITHIUM MINI CHARGER AND LITHIUM DOUBLE CHARGER

### 7.1 POWER SUPPLY

Assemble the power supply as per the instructions in section 7 and connect the lead to the power input port on the charger.



### Charging

No LED	Battery is fully charged.
Flashing LED	Top up charge (Not displayed with NiMH battery)
Solid LED	Battery is charging

The handle can be used at any time during the charging cycle and will automatically resume charging when handle is placed back in the charging well.

When using the Mini charger the handle can be left in place.

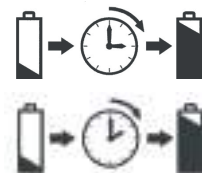


- The instrument must not be used while charging.

## Charging cycle

The Lithium-Ion battery will take approximately 2-3 hours to fully charge. The Lithium-Ion battery will last approximately 2-3 hours on full power.

The NiMH battery will take approximately 1-2 hours to fully charge. The NiMH battery will last approximately 1-2 hours on full power.



## 8. WARRANTY

Your Keeler product is guaranteed for 3 years and will be replaced, or repaired free of charge subject to the following:

- Any fault due to faulty manufacture.
- The instrument and accessories have been used in compliance with these instructions.
- Proof of purchase accompanies any claim.

Please note:

- Batteries are covered by this warranty statement for 1 year only.
- LEDs are covered by this warranty statement for 5 years.
- Bulbs are not covered by this warranty statement.



**The manufacturer declines any and all responsibility and warranty coverage should the instrument be tampered with in any manner or should routine maintenance be omitted or performed in manners not in accordance with these manufacturer's instructions.**

**There are no user serviceable parts in this instrument. Any servicing or repairs should only be carried out by Keeler Ltd. or by suitably trained and authorised distributors. Service manuals will be available to authorised Keeler service centres and Keeler trained service personnel.**

## 9. SPECIFICATIONS AND ELECTRICAL RATINGS

Keeler Direct Instruments and associated power systems are medical electrical instruments. These instruments require special care concerning electromagnetic compatibility (EMC). This section describes the suitability in terms of electromagnetic compatibility of these instruments. When installing or using these instruments, please read carefully and observe what is described here.

Portable or mobile-type radio frequency communication units may have an adverse effect on these instruments, resulting in malfunctioning.

Instrument heads and handles are considered to be inherently EMC benign<sup>1</sup>, with the exception of the GenMed Wall Unit, to which the following table refers, in addition to the Lithium Chargers.

<sup>1</sup> Refer to section 1.4.4 of the Guide for the EMC Directive 2014/30/EU (1st March 2018).

## 9.1 ELECTROMAGNETIC EMISSIONS

### Guidance and manufacturer's declaration – electromagnetic emissions

Keeler Direct Instruments are intended for use in the electromagnetic environment specified below. The customer or the user should assure that they are used in such an environment.

Emissions test		Compliance	Electromagnetic environment – guidance
Chargers and GenMed Wall Unit only	RF emissions CISPR 11	Group 1	Keeler chargers and power systems use RF energy only for their internal function. Therefore, the RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
	RF emissions CISPR 11	Class B	Keeler chargers and power systems are suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes
Harmonic emissions IEC 61000-3-2		Class B	
Voltage fluctuations / flicker emissions IEC 61000-3-3		Complies	

Battery operated Keeler Direct Instruments are considered to be inherently EMC benign<sup>1</sup>, and therefore are not covered by the statements in this section.

<sup>1</sup> Refer to section 1.4.4 of the Guide for the EMC Directive 2014/30/EU (1st March 2018).

## 9.2 ELECTROMAGNETIC IMMUNITY

### Guidance and manufacturer's declaration – electromagnetic immunity

Keeler Direct Instruments are intended for use in the electromagnetic environment specified below. The customer or the user should assure that they are used in such an environment.


Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD). IEC 61000-4-2	± 8 kV contact ± 15 kV air	± 8 kV contact ± 15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst. IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines N/A *± 1 kV for input/output lines	Mains power quality should be that of a typical professional healthcare facility environment.  *GenMed Wall Unit only

Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment – guidance
Surge. IEC 61000-4-5	$\pm 1$ kV line(s) to line(s) $\pm 2$ kV line(s) to earth	$\pm 1$ kV line(s) to line(s) N/A	Mains power quality should be that of a typical professional healthcare facility environment.
Voltage dips, short interruptions and voltage variations on power supply input lines. IEC 61000-4-11	$U_T = 0\%$ 0.5 cycle (0, 45, 90, 135, 180, 225, 270, 315°) $U_T = 0\%$ ; 1 cycle $U_T = 70\%$ ; 25/30 cycles (@ 0°) $U_T = 0\%$ ; 250/300 cycle	$U_T = 0\%$ 0.5 cycle (0, 45, 90, 135, 180, 225, 270, 315°) $U_T = 0\%$ ; 1 cycle $U_T = 70\%$ ; 25/30 cycles (@ 0°) $U_T = 0\%$ ; 250/300 cycle	Mains power quality should be that of a typical professional healthcare facility environment. If the user of the Keeler Direct Instruments requires continued operation during power mains interruptions, it is recommended that the charger be powered from an uninterruptible power supply.
Power frequency (50/60 Hz) magnetic field. IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at a level characteristic of a typical location in a typical professional healthcare facility environment.

Note:  $U_T$  is the a. c. mains voltage prior to application of the test level.

Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment – guidance
			Portable and mobile RF communications equipment should be used no closer to any part of the Digital Camera Assembly, including cables, than the recommended separation distances calculated from the equation applicable to the frequency of the transmitter.
Conducted RF IEC 61000-4-6	6 Vrms 150kHz to 80MHz	6 V	<b>Recommended separation distance</b> $d = 1.2 \sqrt{p}$
Radiated RF IEC 61000-4-3	10 V/m 80MHz to 2.7GHz	10 V/m	$d = 1.2 \sqrt{p}$ 80MHz to 800 MHz $d = 2.3 \sqrt{p}$ 800MHz to 2.7GHz



Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment – guidance
			<p>Where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey<sup>1</sup>, should be less than the compliance level in each frequency range.<sup>2</sup></p> <p> Interference may occur in the vicinity of equipment marked with this symbol.</p>

Note 1: At 80MHz and 800MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

1 Field strengths from fixed transmitters, such as base stations (cellular / cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Direct Instruments are used exceeds the applicable RF compliance level above, the Direct Instruments should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orientating or relocating the Direct Instruments.

2 Over the frequency range 150kHz to 80 MHz, field strengths should be less than 10 V/m.

### 9.3 RECOMMENDED SAFE DISTANCES

#### Recommended separation distances between portable and mobile RF communications equipment and Keeler Direct Instruments.

Keeler Direct Instruments are intended for the use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of Keeler Direct Instruments can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Keeler instruments as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80MHz $d = 1.2\sqrt{p}$	80MHz to 800MHz $d = 1.2\sqrt{p}$	800MHz to 2.7GHz $d = 2.3\sqrt{p}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance in metres (m) can be determined using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80MHz and 800MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

## 10. TECHNICAL SPECIFICATIONS

The Ophthalmoscope / Retinoscope / Otopscope, the power supply (EP29-32777) with its charging dock (1941-P-5289 and 1941-P-5326) together constitute a Medical Electrical System as defined in EN/IEC 60601-1.


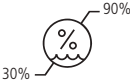
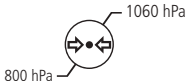
### Power Supply

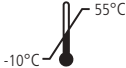
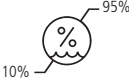
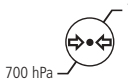

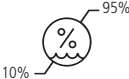
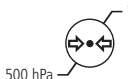
<b>Input mains data</b>	100-240V – 50/60Hz
<b>Power supply rating</b>	12V: 2.5amps
<b>Operation</b>	Maximum 15 minutes ON Minimum 10 minutes OFF
<b>Classification:</b>	Class II equipment Type B protection against shock

### Instrument heads and handles

<b>Input voltage (DC)</b>	3V 2xAA Alkaline Batteries - BLUE 3.75V Lithium-Ion Rechargeable Battery - RED (EP39-18918) 3.65V NiMH Rechargeable Battery - Black (1919-P-7149)
---------------------------	---

### Environmental Conditions:

USE	
	
	
Shock (without packing)	10 g, duration 6 ms

STORAGE CONDITIONS	
	 
TRANSPORT CONDITIONS	
	 
Vibration, sinusoidal	10 Hz to 500 Hz: 0.5g
Shock	30 g, duration 6 ms
Bump	10 g, duration 6 ms

## 11. ACCESSORIES AND SPARES

Item	Part Number
Spec/Vista halogen bulb 3.6V (Pack of 2)	1011-P-7034
Spec/Vista halogen bulb 2.8V (Pack of 2)	1011-P-7042
Spec/Vista LED bulb 2.8V/3.6V (Pack of 1)	1011-P-7229
<b>Ophthalmoscopes Standard</b>	
Std Oscope Halogen Bulb 2.8V (Pack of 2)	1015-P-7031
Std Oscope Halogen Bulb 3.6V (Pack of 2)	1015-P-7023
Std/Pract/Prof Xenon Bulb 2.8V (Pack of 2)	1011-P-7106
Std/Pract/Prof Xenon Bulb 3.6V (Pack of 2)	1011-P-7114
Ophth LED Assy	1011-P-5610
<b>Ophthalmoscopes Practitioner</b>	
Fo Oscope Halogen Bulb 2.8V (Pack of 2)	1015-P-7066
Fo Oscope Halogen Bulb 3.6V (Pack of 2)	1015-P-7058
Std/Pract/Prof Xenon Bulb 2.8V (Pack of 2)	1011-P-7106
Std/Pract/Prof Xenon Bulb 3.6V (Pack of 2)	1011-P-7114
Ophth Led Assy	1011-P-5610
<b>Fibre-Optic Oscope</b>	
FO Oscope Halogen Bulb 2.8V (Pack of 2)	1015-P-7066
FO Oscope Halogen Bulb 3.6V (Pack of 2)	1015-P-7058

Item	Part Number
Std/Pract/Prof Xenon Bulb 2.8V (Pack of 2)	1011-P-7106
Std/Pract/Prof Xenon Bulb 3.6V (Pack of 2)	1011-P-7114
1011-P-5610 Ophth LED Assy	
<b>Pocket</b>	
Std Otscope Halogen Bulb 2.8V (Pack of 2)	1015-P-7031
Pocket Ophth Halogen Bulb 2.8V (Pack of 2)	1011-P-7050
<b>Other – Chargers</b>	
Lithium Double Charger	1941-P-1368
Lithium Mini Charger	1941-P-1341
3.6V Lithium Battery	EP39-18918
<b>Other – Colour coded grips</b>	
Slimline Handle Sleeve – Pink	1901-P-7028
Slimline Handle Sleeve – Green	1901-P-7036
Slimline Handle Sleeve – Blue	1901-P-7044
Slimline Handle Sleeve – Black	EP29-05365
Slimline Handle Sleeve Assorted Colours	1901-P-7052
<b>Other – Specula – Jazz Ultra</b>	
Jazz 2mm Reusable Specula (Pack of 10)	1514-P-7036
Jazz 2.5mm Reusable Specula (Pack of 10)	1514-P-7044
Jazz 3mm Reusable Specula (Pack of 10)	1514-P-7052
Jazz 4mm Reusable Specula (Pack of 10)	1514-P-7060
Jazz 5mm Reusable Specula (Pack of 10)	1514-P-7079
Jazz 2mm Specula (Pack of 100)	1514-P-7087
Jazz 2.5mm Specula (Pack of 100)	1514-P-7095
Jazz 3mm Specula (Pack of 100)	1514-P-7108
Jazz 4mm Specula (Pack of 100)	1514-P-7116
Jazz 5mm Specula (Pack of 100)	1514-P-7124

## 12. PACKAGING AND DISPOSAL INFORMATION

### Disposal of old electrical and electronic equipment



This symbol on the product or on its packaging and instructions indicates that this product shall not be treated as household waste.

To reduce the environmental impact of WEEE (Waste Electrical Electronic Equipment) and minimise the volume of WEEE entering landfills we encourage at product end of life that this equipment is recycled and reused.

**If you need more information on the collection reuse and recycling then please contact B2B Compliance on 01691 676124 (+44 1691 676124). (UK only).**

**Any serious incident that has occurred in relation to the device must be reported to the manufacturer and the competent authority of your Member State.**

## Contact



### Manufacturer

Keeler Limited  
Clewer Hill Road  
Windsor  
Berkshire  
SL4 4AA UK

**Freephone** 0800 521251

**Tel** +44 (0) 1753 857177

**Fax** +44 (0) 1753 827145

### USA Sales Office

Keeler USA  
3222 Phoenixville Pike  
Building #50  
Malvern, PA 19355 USA  
**Toll Free** 1 800 523 5620  
**Tel** 1 610 353 4350  
**Fax** 1 610 353 7814

### India Office

Keeler India  
Halma India Pvt. Ltd.  
Plot No. A0147, Road No. 24  
Wagle Industrial Estate  
Thane West – 400604, Maharashtra  
INDIA  
**Tel** +91 22 4124 8001

### China Office

Halma China Group  
名称: 沃迈(上海)机电有限公司  
地址: 上海市闵行区金都路1165弄  
123号23幢一号厂房三层B座  
电话: 021-6151 9025



Visiometrics, S. L., Vinyals, 131  
08221 Terrassa, Spain

EP59-11234 Issue 9

Date of Issue 12/05/2021



**Keeler**  
– A world without vision loss –